




- Buckle system easily accommodates highly compressible thigh tissue
- Light weight, strength tested, durable materials
- Customizable socket to maximize comfort for each individual patient
- Total contact giving a more stable gait and better range of motion
- Unique subischial design provides better sitting comfort
- Wide range of adjustments that can be made by the user
- Patients can adjust their fit any time without using socks
- A specially designed light weight offset offers rotation and various degrees of flexion of the socket to facilitate the alignment process

Specifications

	Mobility Grade	K0 - K3
	Side	Left (L), Right (R)
	Distal Circumference (cm)	25 - 44
	Limb Length (cm)	15 - 30 +
	Product weight (Includes lock and pyramid)	3.26 lbs/ 1.48 kg
	Weight limit	136 kg / 300 lbs
	Buckles	Right or Left
	Closure System	Cable or Spool
	Water Capability	Waterproof
	Warranty	1 year
	Device Models	TF300R / TF300L

Prosthesis Kit

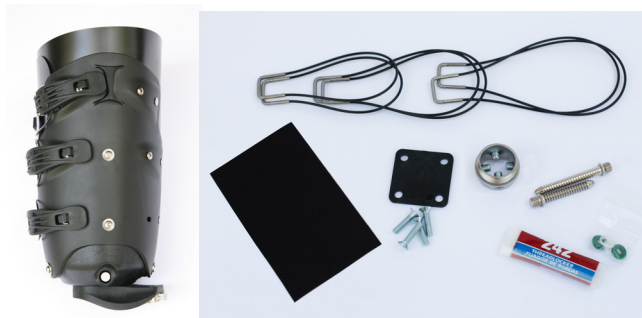

- Adjustable socket with flexible inner liner
- Rotating offset with attached rotating pyramid receiver
- Non-slip grip material
- Five cables (45cm-25cm)
- Two locking pins
- 2 degree wedge with longer screws
- Loctite 242
- Loop velcro
- Utility cord
- Three spool conversion kits (spool, cord and screws)

Call Joshua Mullins at 301-821-5414 or email jmullins@ifitprosthetics.com for more information.

iFITProsthetics.com | Follow us:

iFIT Prosthetics is a Veteran-Owned Small Business (VOSB)

Patent #8470050, #8491667, #8845755, #10398577, #10806608, #11382775, #11766343. ©2024 iFIT Prosthetics, LLC. All rights reserved. Products mentioned herein are trademarks or registered trademarks of iFIT Prosthetics, LLC.



The Right Fit, Right Now™